BCH 462- Biotechnology & Genetic engineering [Practical] Lab (4) Western Blot


Immunoassay

- Antigens [Ag]:
- A substance that when introduced into the body stimulates the production of an **antibody**.
- Antigens include toxins, bacteria, foreign blood cells, and the cells of transplanted organs.

Pause and Think how to prevent the body from rejecting an organ transplant?


- Antibody [Ab]:
- Antibodies are large Y-shaped glycoproteins.
- They are produced by the immune system to identify and neutralize foreign objects (antigens).
- Immunoassay:

A test that uses highly specific and selective <u>antigen-antibody reactions</u> forming antibody and antigen complexes [**immuno-complexes**] <u>as a means of generating measurable results.</u>

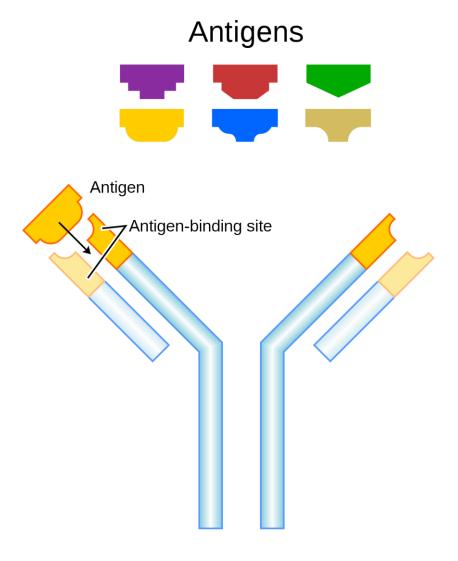

Antibody

Antibody and antigen complexes [immuno-complexes]

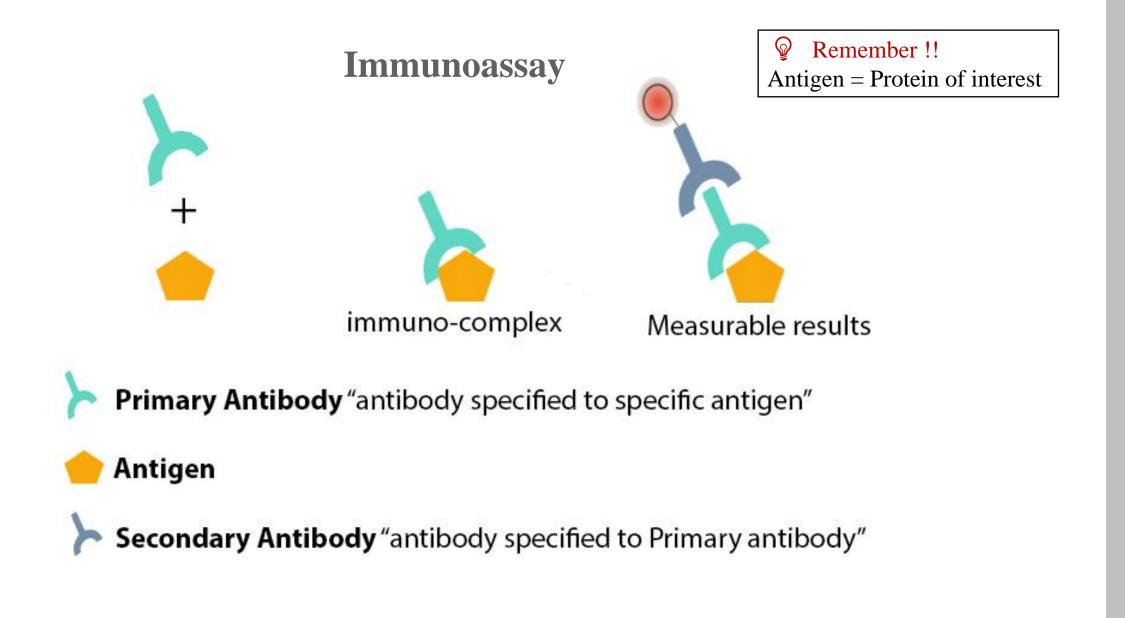


Blotting

 Blotting is the <u>transfer of macromolecules</u> such as nucleic acids and proteins to <u>solid-phase</u> <u>membranous support</u> (for instance, a nitrocellulose, PVDF or nylon membrane).

• Types:

- 1. Southern blotting used to detect DNA.
- 2. Northern blotting used to detect RNA.
- 3. Eastern blotting; used to analyze protein post translational modifications (PTM).
- 4. Western blotting used to detect Proteins.

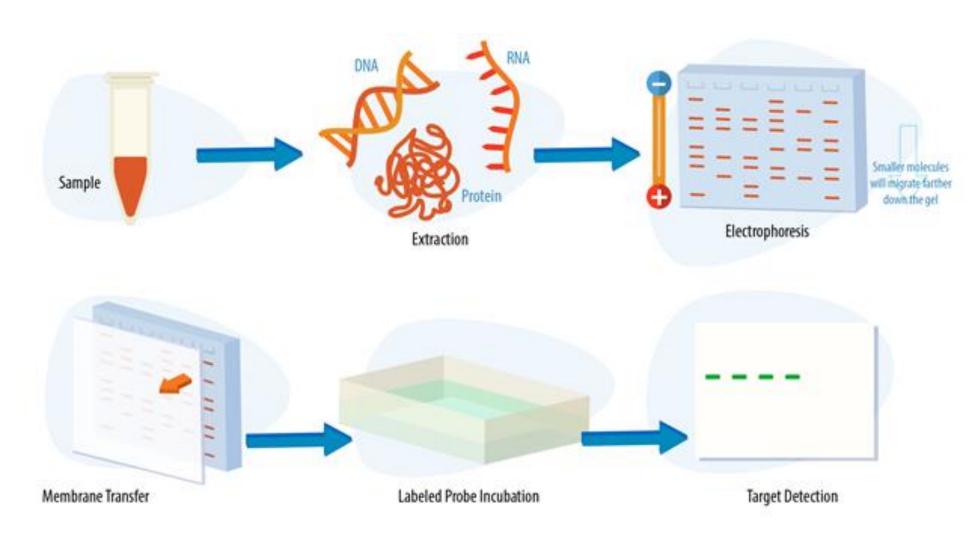

Overview of blotting technique

Western blot

- Also called **protein immunoblot**.
- It is a widely used immunoassay technique.
- Used to identify specific proteins [antigens] in a sample of tissue homogenate or extract, based on their ability to bind to [antibodies] resulting in color/florescence indicate the presence of this specific protein.
- Applications:
- 1. Analyzing, identifying target proteins and estimating their molecular weight.
- 2. To **compare** the amounts of a protein of interest among different samples.
- 3. Used in clinical laboratories for assisting **identification** of certain antigen proteins (pathogen or biomarker).
- 4. Used to detect **changes in protein expression** under different biological conditions (e.g. in disease, stress, etc.).

Practical Part 🧬

Practical part

- Aims:
 - To understand how proteins (antigens) can be analyzed using antibodies raised against these proteins by <u>immunoblotting technique</u>.
 - To perform the steps of western-blot technique to detect <u>tubulin protein</u>.

Practical part

- The mixture of proteins is **separated** based on **molecular weight**.
- These results are then <u>electro-transferred</u> to solid support producing a band for each protein.
- The transferred protein is detected by incubating the membrane with specific primary antibody to the protein
 of interest, secondary antibody labelled with an enzyme or fluorophore which target the primary antibody,
 and substrate which in the end you will get colored product.
- The **color/florescence** indicates <u>the presence</u> of the protein of interest.
- The thickness of the band corresponds to the amount of protein present.
- Thus, the molecular weight and amount of the desired protein can be characterized from a complex mixture of proteins by western blotting.

Western blot performing steps


The technique uses three elements to accomplish this task




1. 1st phase (SDS-PAGE):

Separation sample mixture by size using SDS-PAGE.

2. 2nd phase (Electro-blotting):

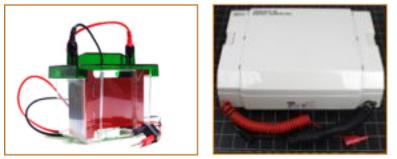

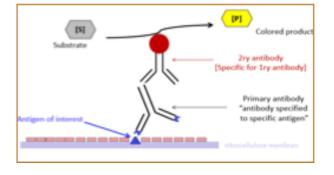


Transfer to a solid support (**electro-blotting**) by transferring the proteins bands from the gel to the membrane.


3. 3rd phase (Marking target protein to visualize):
Marking target protein using a proper primary and secondary antibody to visualize.

1st phase (SDS-PAGE)


• A protein sample is subjected to polyacrylamide gel electrophoresis.

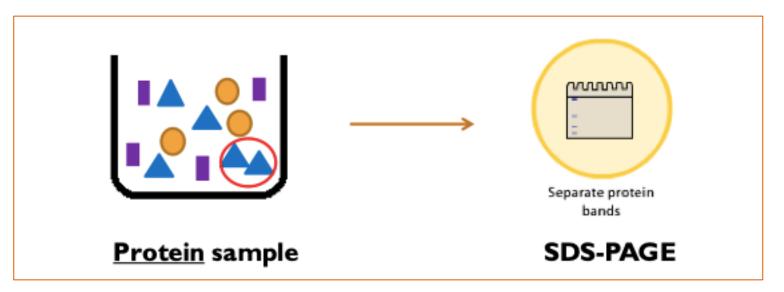

- To confirm the separation of samples (since separated proteins are colorless) use:
 - 1. Replica of the gel and stain it as usual [with Coomassie brilliant blue R-250].
 - 2. Using a pre-stained marker.
 - 3. Reversible staining by Ponceau S.

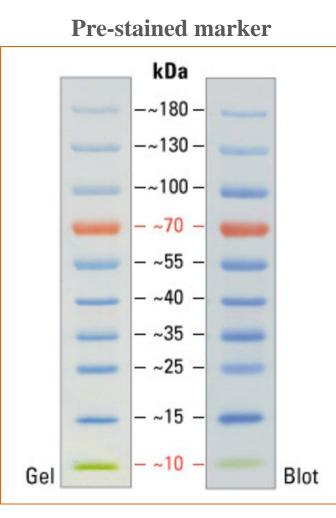


Pre-stained Protein Ladder is a mixture of 10 blue-, orange-, and green-stained proteins (10 to 180 kDa) for use as size standards in protein electrophoresis (SDS-PAGE) and western blotting.

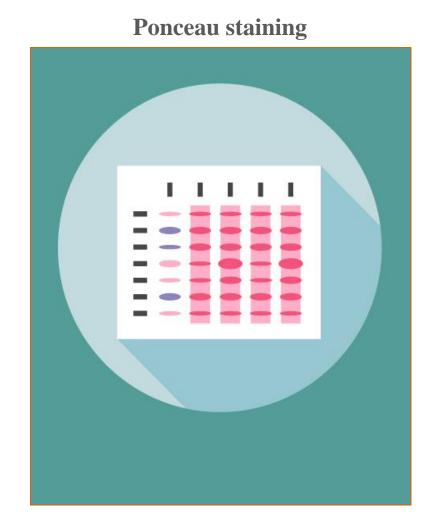

Ponceau staining is a washable light red colored dye, that may be used to prepare a stain for rapid detection of protein bands on nitrocellulose or polyvinylidene fluoride (PVDF) membranes (Western blotting).

2nd phase (Electro-blotting)


 After that the gel is placed over a sheet of PVDF, the protein in the gel is electrophoretically transferred to the Polyvinylidene fluoride (PVDF) membrane. "transfer step [Electroblotting]".

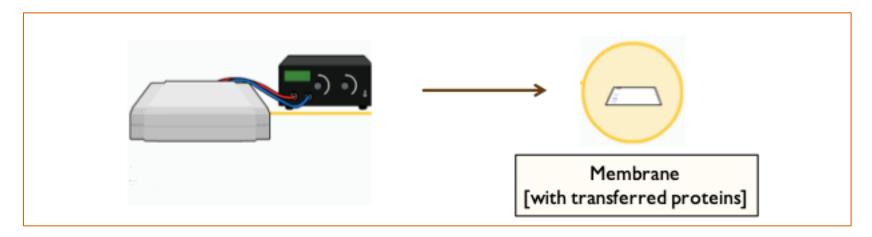

Methods of transfer:

- 1. Wet method \rightarrow Most common transfer method, best for proteins >100kDa.
- 2. Semi-wet \rightarrow Quick but less efficient.
- 3. **Dry** \rightarrow Quicker, no transfer buffer required.


Transfer sandwich

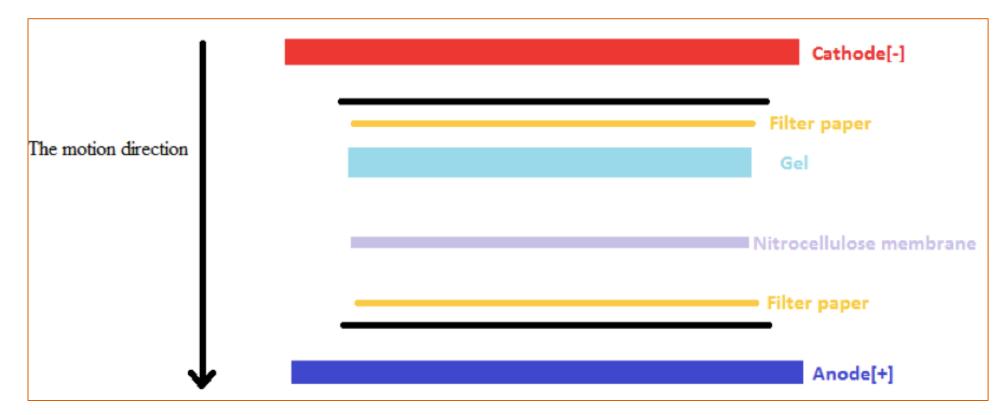


- Because the samples in the gel are [-ve] charged, the applied electric current will facilitate their transferring to nitrocellulose membrane, the samples will move toward the Anode[+].
- Also the capillary action has its effect in the movement of the samples from the gel to the nitrocellulose membrane.
- Note that: [the filter papers, gel and PVDF membrane will be soaked in transfer buffer].

3rd phase (Marking target protein to visualize)

• The PVDF is then soaked in <u>blocking buffer</u>.

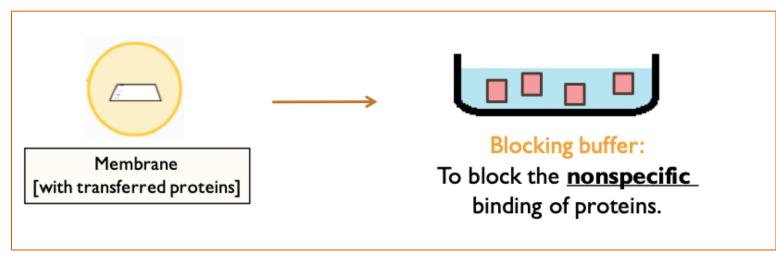


• The PVDF is then incubated with the <u>specific primary antibody</u> for the protein of interest.

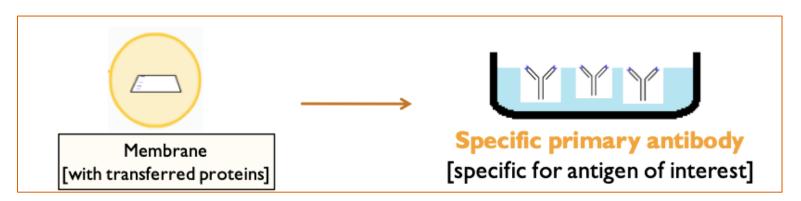

3rd phase (Marking target protein to visualize)

- The PVDF is then washed and incubated with a <u>secondary antibody</u>, which is specific for the first antibody [primary_antibody].
- 2^{ry} AB is conjugated with enzymes or fluorescent that give a subsequent reaction with an applied reagent, leading to a coloring or emission of light, enabling detection.

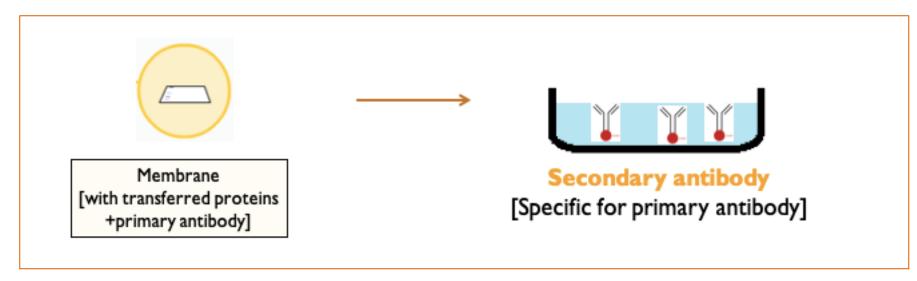

• The color produced indicate the **presence** of the **antibody-antigen [Ab-Ag]** complex.

3rd phase (Marking target protein to visualize)

Types of western blot detection methods "detection step":

1- Colorimetric

Secondary antibody labeled with **enzyme** conjugates, such as alkaline phosphatase (AP) or horseradish peroxidase (HRP), when provided with a chromogenic substrate, will cause a color reaction.


2- Chemiluminescence


The **enzyme** attached to the secondary antibody triggers a reaction with <u>a luminescent substrate that produces</u> <u>light as a by-product.</u>

3- Fluorescence

Uses secondary antibodies that are conjugated to <u>specific fluorophores</u> that can absorb and emit light within a <u>range of wavelengths</u>, so no additional substrate is necessary.


Western blot detection methods:

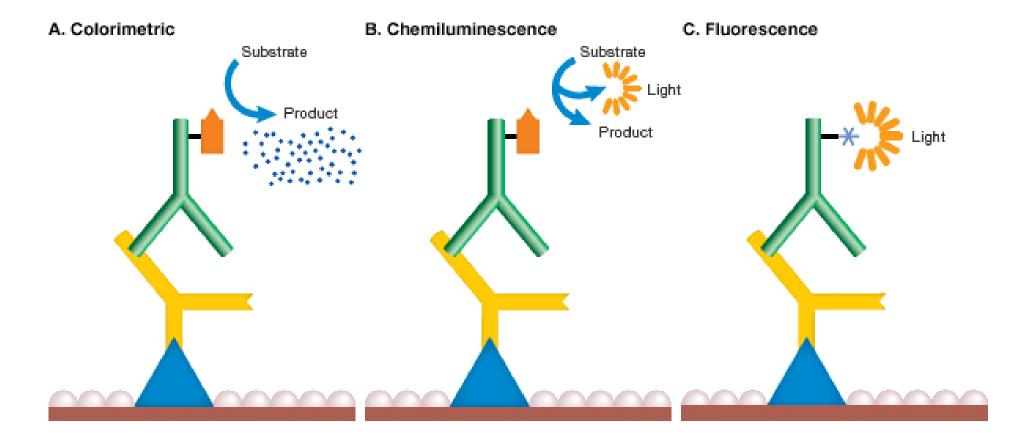

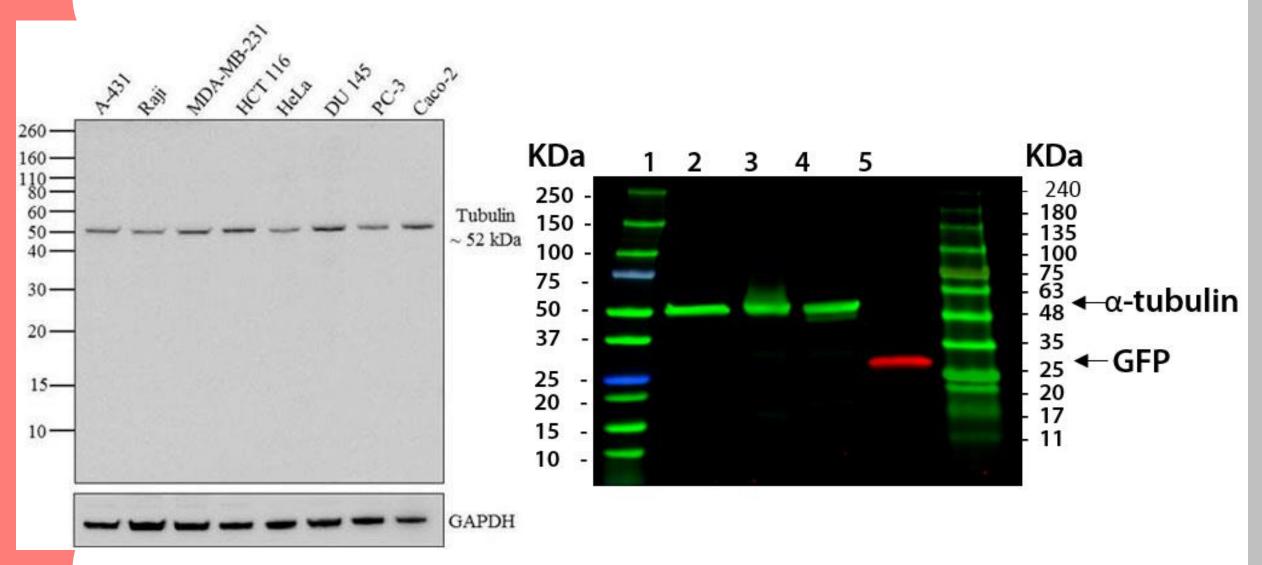

- Thus the molecular weight and amount of the desired protein can be characterized from a complex mixture of proteins by western blotting.
 - Pause and Think What is the Mwt of tubulin? Which cell-line has the highest tubulin concentration? What types of test is Western blot? quantitative OR qualitative?

Detection of specific protein using Western blot

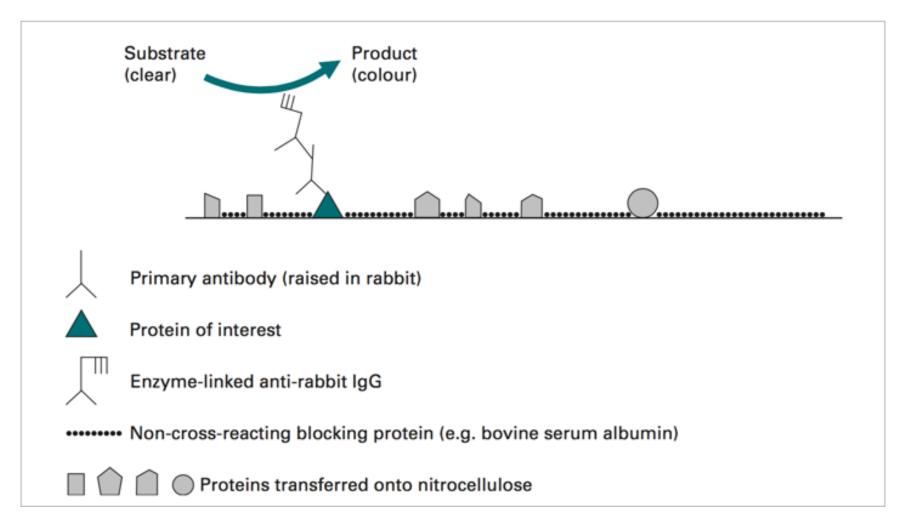

Supporting materials:

Performing western blot:

http://www.youtube.com/watch?v=VgAuZ6dBOfs